Wisd Instruments



Magnetic Stirrer with Hotplate MSH-500D



Magnetic Stirrer with Hotplate Highest performance - highest convenience - highest user safety!

Our products are well known in almost 100 countries around the world for their performance, outstanding quality and highest user safety at a very competitive price. The new Digital Magnetic Stirrer with Hotplate was designed to fullfill the highest requirements in the laboratory.



Highest performance

- powerful: stirring capacity up to 20 l, temperature up to 500 °C, stirring rate up to 1.500 rpm
- protection against corrosion due to the smart design
- -superior heat and chemical resistant full ceramic plate
- -strong but quiet motor, made for a long life time
- precise temperature by digital feedback control

Highest convenience


- intuitive and very easy-to-use control elements
- illuminated digital display for easy setting
- •easy-to-clean full ceramic plate



Highest user safety

- hot top indicator for temperatures over 50°C
- high quality insulation to protect user and device against heat damage
- two independent internal sensors for accurate temperature and overheating protection



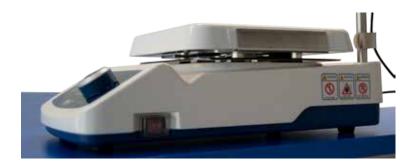

500°C heating power

Protection

Two features protect the device from internal corrosion and prolong its life time: the **closed coverage type** and the **intelligent fluid management system**, which leads off accidentally spilled liquids away from the interior of the device.

Surface

The superior full ceramic plate achieves up to 500°C, is resistant against aggressive chemicals and easy-to clean.

Housing

The **die-cast case** with **high-quality insulation** ensures **user safety** during operation at high temperatures. The **minimized heat transfer** from the plate to the inside also ensures a **longer life time**.

(S witeg

Accessories

The external temperature sensor allows accurate temperature due to digital feedback control system.

Motor

The brushless shade motor ensures strong and precise stirring, minimized noise and a long life time. Due to the high-performance magnet, up to 20 I of water can be stirred easily.

Safety first

Hot top indicator lamp illuminates when temperature on the plate is over 50°C.

Temperature control & safety

Two independent internal temperature sensors assure **accurate temperature control** and **prevent from malfunction** or danger caused by overheating.

Controller

New jog-dial "turn & push" and touch button are designed for intuitive, easy and direct use.

Wisd Instruments



Magnetic Stirrer with Hotplate MSH-500D

Your advantages at a glance:

- highest quality at a very competitive price
- powerful: up to 20 I stirring capacity, up to 500 °C, up to 1 500 rpm
- highest user safety to prevent accidental burning
- dual protection against internal corrosion
- precise temperature control by digital feedback
- user-friendly design for highest convenience
- high-quality materials and stable design for a long life time
- CE label
- 24 months warranty

Order information



 Device
 Order number

 Magnetic stirrer with hotplate MSH-500D
 DH.WMH05020


 Magnetic stirrer with hotplate MSH-500D Set including:
 DH.WMH04020

 • temperature probe
 DH.WMH04020

 • stand rod
 DH.WMH04020

 • clamp and clip for temperature probe
 DH.WMH04020

Magnetic stirrer MS-500D	DH.WMS04022
Hotplate HP-500D	DH.WHP05020
Hotplate HP-500D Set including: • temperature probe • stand rod	DH.WHP04020

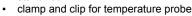



MS-500D

witeg Labortechnik GmbH Am Bildacker 16 97877 Wertheim Germany

Tel.: +49 (0)9342 / 9301-0 Fax: +49 (0)9342 / 9301-77 www.witeg.de info@witeg.de