


The most comfortable way for all pipetting applications

Looking for the most comfortable way of pipetting? The electronic microliter pipette **WITOPET EP** is the right product for you. The motorized pipetting process has at least two major advantages compared with manual pipetting: it avoids pipetting errors user induced and it prevents the risk of a cramp after prolonged pipetting phases. **WITOPET EP** is available as single- and multi-channel pipette.

Digital display

The user-friendly display indicates set volume, set function, pipetting step, battery status etc.

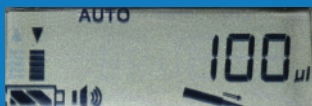

Volume and function setting

The buttons for volume and function setting are color-coded to avoid confusion.



Push-button


Draw up and release of the set volume occurs with minimal effort for highest user comfort.

Tip ejector

The tip ejector releases the used tip easily.

Pipetting unit

The pipetting unit can be turned around (360°) for increased flexibility.

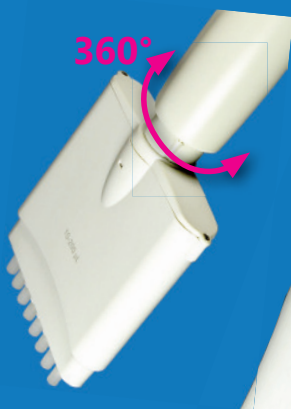

Body

A high precision stepping motor ensures absolute correct automatic pipetting without any mistake.

Tip cone

The tip cone is suitable for all established brands of pipette tips.

Pipette support

suitable for **WITOPET EP** pipettes
5 406 105





Electric Pipette

Your advantage at a glance:

- ideal for prolonged pipetting
- minimal pipetting forces required
- automatic pipetting process for highest comfort and accuracy
- multiple functions (auto-pipetting, mixing, multiple dispensing, sequential aspirating / dispensing)
- 9 memory settings
- 5 speed settings for different viscosities
- user friendly digital display
- 360° rotatable pipetting unit for more flexibility (multi-channel only)
- including Li-ion 3,6V 1000mAh battery
- serial number and CE label
- individual certificate according to DIN EN ISO 8655-2
- DE-M marked according to the German calibration law

WITOPET EP digital single- and multi-channel pipette

	Volume	Channels	Division	Volume tested	Systematic errors	Random errors	Order number
	μl		μl	μl	$\pm \mu\text{l}$	$\leq \mu\text{l}$	
single-channel	2 - 20	1	0.1	2 10 20	0.2	0.1	5 402 302 S
	10 - 200	1	1.0	20 100 200	1.6	0.6	5 402 320 S
	100 - 1000	1	10.0	200 500 1000	8	3	5 402 390 S
multi-channel	2 - 20	8	0.1	2 10 20	0.4	0.2	5 402 402 S
	10 - 200	8	1.0	20 100 200	3.2	1.2	5 402 420 S
Pipette support for WITOPET EP							5 406 105